

FAUCET-MOUNTED EYEWASH KIT



Application: Faucet-mounted eyewashes convert any faucet into an emergency eyewash station without interfering with normal faucet operation. Faucet-mounted eyewash is ideal for use with gooseneck faucets, and can be installed at any sink, close to where accidents might occur. In an emergency, the unit is quickly located and activated, and provides an unlimited supply of potable water for rinsing the user's eyes.

Outlet Heads: Outlet heads are mounted 5" apart and deliver a soft, aerated flow of water. Heads angle forward and inward toward user. Angle of heads is adjustable to permit full coverage and avoid splashing. Furnished with float-off dust covers to protect outlet heads.

Valve: Forged brass diverter valve. Pull knob to activate eyewash; water pressure holds eyewash in operation, leaving user's hands free. Push knob or turn off faucet to return to normal faucet operation.

Inlet: Body has 55/64"-27 female thread. Furnished with three adaptors (15/16"-27, 13/16"-27 and 3/8" NPS) for installing on most commonly used faucets, including laboratory-type faucets.


Outlet: Furnished with removable aerator on bottom.

Quality Assurance: Each unit is completely assembled and water tested prior to shipment.

IMPORTANT: Faucet-mounted eyewashes require two motions to operate (turn on water, pull knob to activate eyewash flow).

Faucet-mounted eyewashes should be used with cold or warm water only. Use of hot water might cause scalding..


FAUCET-MOUNTED EYEWASH KIT

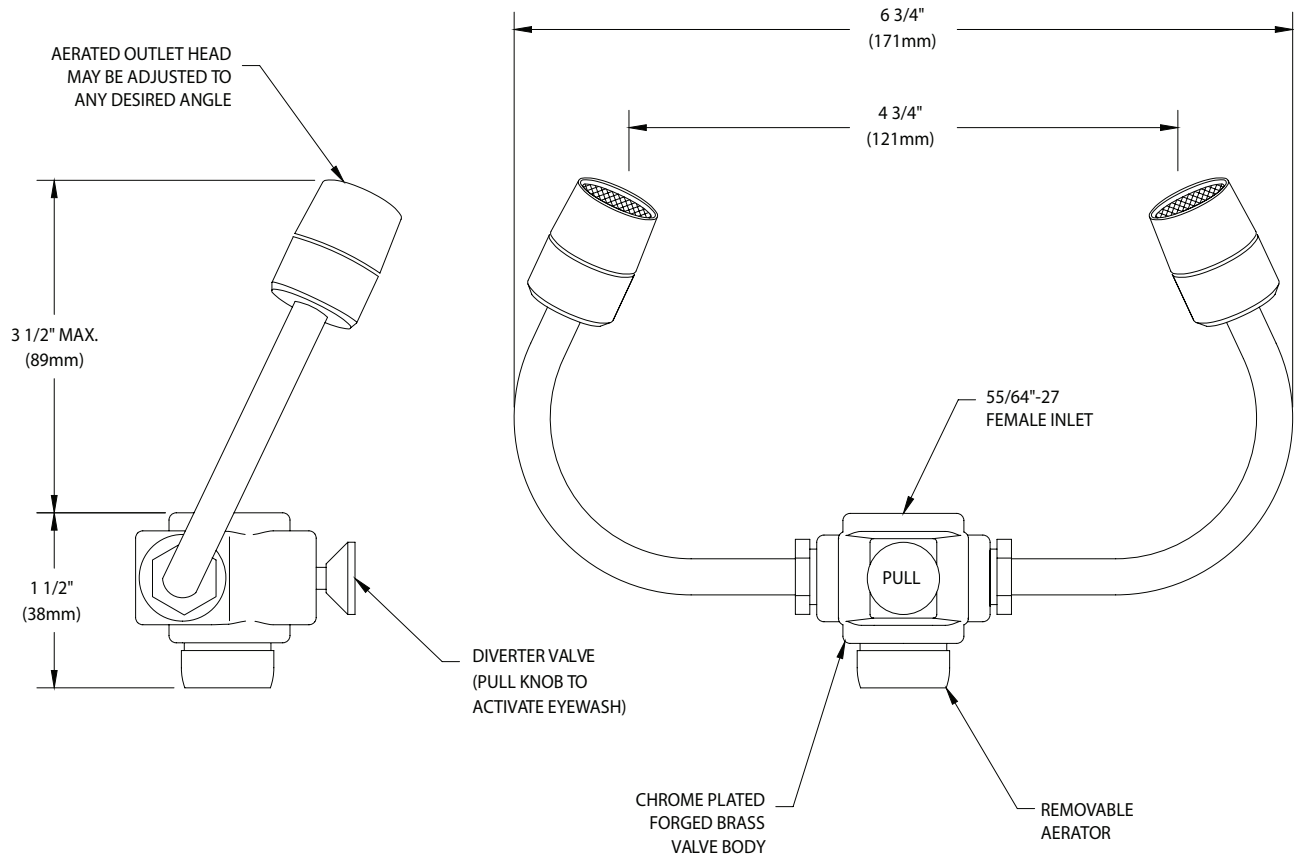

15/16"-27 ADAPTOR



13/16"-27 ADAPTOR

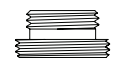

3/8" NPS ADAPTOR

NOTES:

1. USE COLD OR WARM WATER ONLY WITH FAUCET-MOUNTED EYEWASH.
2. THIS UNIT IS INTENDED AS A SUPPLEMENT TO, BUT NOT A REPLACEMENT FOR, DEDICATED PLUMBED-IN EYEWASH EQUIPMENT.