

TECLAB IS RECOGNIZED AS A LEADER IN THE DESIGN, MANUFACTURE AND INSTALLATION OF LABORATORY FURNITURE, ERGONOMIC WORKSTATIONS, WORKBENCH SYSTEMS, AND STORAGE EQUIPMENT.

Whether you're looking to create a single workstation or an entire laboratory facility, Teclab can supply all the essential components for a highly efficient and productive work space.

We are an **American** company that's been in the laboratory furniture business since **1982**, supplying top companies in the energy, pharmaceutical, healthcare, biotech and research fields as well as government and educational facilities. These organizations rely on Teclab to supply the best products and workmanship.


Customer service has always been a top priority at Teclab and our extremely knowledgeable staff can also assist you in planning your laboratory. From room layout to installation, our in-house designers and engineers use their decades of laboratory design experience to customize your application needs.

You can call us at **269.372.6000** and speak one on one with a helpful and friendly member of Teclab's planning and design team or send us an email at **info@teclab.com** and start creating your new laboratory.

**HIGHEST QUALITY. BEST BUILT.
MOST DURABLE. MOST VERSATILE.**

The extraordinary fit and finish of every Teclab product reflects that devotion to providing what customers want and need. Customers see our quality in clean, functional lines. They hear it in the sound of a drawer sliding effortlessly closed. They feel it in the solidity of a well made unit. And they experience it in the attention to the smallest of details, our design assistance, installation and after-sale service.

- **MANUFACTURED IN-HOUSE USING FURNITURE GRADE COLD ROLLED STEEL.**
- **10 YEAR WARRANTY ON ALL TECLAB BUILT PRODUCTS***
- **MODULAR – BUILD YOUR OWN LAYOUT**
- **BUILT TO ORDER FOR EASY CUSTOMIZATION**

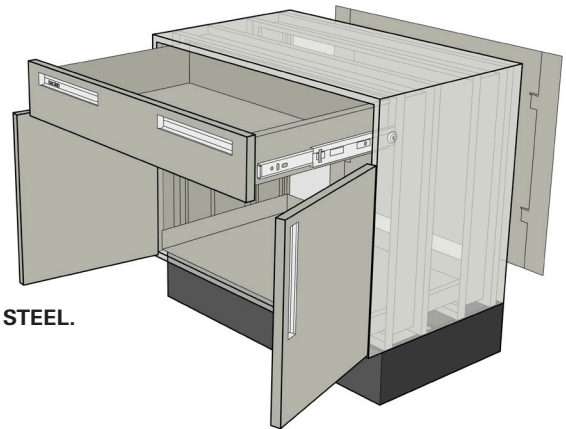

DRAWERS:

Smoothly opening drawers put contents at your fingertips without sagging. Precision-made drawer slides, consisting of an outer channel, intermediate channel and center track move freely on a continuous series of ball bearings for effortless operation. Accurately machined grooves assure precise fit to eliminate any side-play between members. Special rollers facilitate the disengagement and re-entry of the integral chassis and track section into the remaining side members and subsequently into the cabinet. Heavy-duty cold rolled steel drawer construction combined with the precision slide mechanism results in related load capacities far in excess of requirements.

HANDLES:

Teclab inset handles (shown above) are extruded aluminum with an anodized finish (11" w x 1-1/4" h x 1" d). Wetlab external handles are stainless steel (5" w x 1-3/8" d x 3/8" dia.).

HINGES:

Concealed soft-close hinges. Tension closed (no latches required).

LOCKS:

Security key-locks can be installed in any Teclab drawer cabinet at the factory, or at any time with knockouts provided.

DOORS:

Heavy-duty doors feature double wall construction of heavy duty cold rolled steel with front panel and back liner for sound-proofing. Firm, quiet closing is assured by the use of rubber bumpers installed in every door liner.

INTERNAL FRAMES:

High quality Base Cabinets with rugged twin internal frames form the foundation for Teclab Work Station superiority.

REMOVABLE BACK PANELS:

All Teclab base cabinets feature a removable back panel. External on technical style cabinets (shown above) and internal on wetlab style cabinets.

*Parts not manufactured by Teclab, such as height adjustment mechanisms, casters, upholstery, etc. have their own warranties.



TAU - v.1.1