


BASE CABINET CONSTRUCTION DETAILS

HIGHEST QUALITY. BEST BUILT. MOST DURABLE. MOST VERSATILE.


Take your pick. They're all descriptions that have become synonymous with Teclab in our 42 year History. But our dedicated commitment to building the best is more than just talk.

The extraordinary fit and finish of every Teclab product reflects that devotion to providing what customers want and need. Customers see our quality in clean, functional lines. They hear it in the sound of a drawer sliding effortlessly closed. They feel it in the solidity of a well made unit. And they experience it in the attention to the smallest of details, our design assistance, installation and after-sale service.

MANUFACTURED IN-HOUSE. MODULAR — BUILD YOUR OWN LAYOUT BUILT TO ORDER FOR EASY CUSTOMIZATION

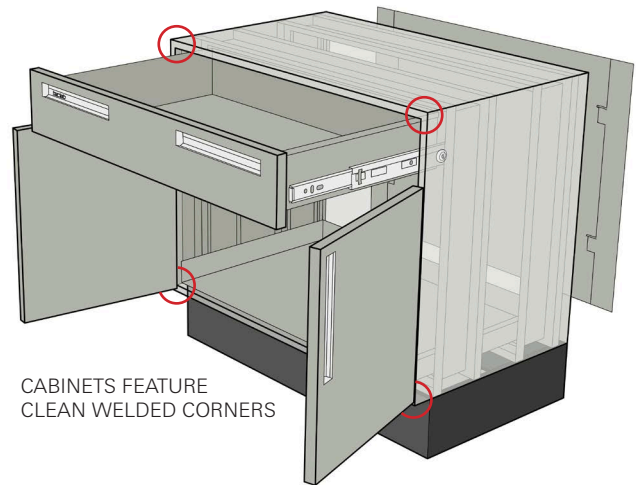

CABINETS FEATURE
CLEAN WELDED CORNERS

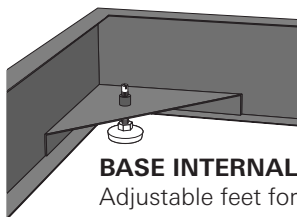

BASE INTERNAL DETAIL:
Adjustable feet for uneven floors.

STEEL GAUGES:

- 20 Ga.** CABINET SHELLS, DRAWERS, DOOR FACES
- 18 Ga.** BASES, INTERNAL FRAMES, DOORS (INTERNAL)
- 12 Ga.** BASE SUPPORTS

Also applies to our Stainless Steel Cabinets.

DRAWERS:

Smoothly opening drawers put contents at your fingertips without sagging. Precision-made drawer slides, consisting of an outer channel, intermediate channel and center track move freely on a continuous series of ball bearings for effortless operation. Accurately machined grooves assure precise fit to eliminate any side-play between members. Special rollers facilitate the disengagement and re-entry of the integral chassis and track section into the remaining side members and subsequently into the cabinet. Heavy-duty cold rolled steel drawer construction combined with the precision slide mechanism results in related load capacities far in excess of requirements.

HANDLES:

Teclab inset handles (shown above) are extruded aluminum with an anodized finish (11" w x 1-1/4" h x 1" d). Wetlab external handles are stainless steel (5" w x 1-3/8" d x 3/8" dia.).

HINGES:

Concealed soft-close hinges. Tension closed (no latches required).

LOCKS:

Security key-locks can be installed in any Teclab drawer cabinet at the factory, or at any time with knockouts provided.

DOORS:

Heavy-duty doors feature double wall construction of heavy duty cold rolled steel with front panel and back liner for sound-proofing. Firm, quiet closing is assured by the use of rubber bumpers installed in every door liner.

INTERNAL FRAMES:

High quality Base Cabinets with rugged twin internal frames form the foundation for Teclab Work Station superiority.

REMOVABLE BACK PANEL:

All Teclab base cabinets feature a removable back panel. External on technical style cabinets (shown above) and internal on wetlab style cabinets.