


**BENCH HIGH BASE CABINET C-1504**

Teclab Base Cabinets are constructed of only top quality, heavy-duty, furniture-grade, cold-rolled steel, double-braced and welded at all critical joints; twin full-perimeter frames formed of heavy gauge steel, welded to provide strength and rigidity; heavy-duty drawers that glide effortlessly on precision steel ball bearings; cabinet doors that open and close with authority and precision using self-closing hinges; brushed aluminum pulls that are recessed for clean lines and unobtrusive operation; durable, textured paint finish, in a wide selection of mar-resistant, easy-to-clean colors.

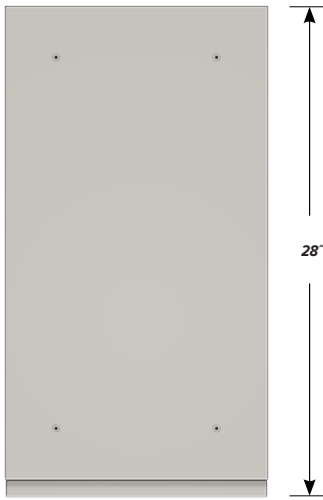

**TOP**

**C-1504 FEATURES**

**One "C" drawer, one "A" door cabinet with LH pull and sliding bottom tray.**

Cabinet Dimensions: 15" wide, 33 3/4" high, 28" deep

**"C" Drawer**

Drawer Face Dimensions: 15" wide, 13 1/2" high, 1" thick.  
 Drawer Body Inside Dims.: 12" wide, 9 1/2" high, 25 1/2" deep.  
 Weight Capacity: 75 lbs.


*All 15" wide "C" drawers accept standard hanging files.*

**"A" Door with Left Hand Pull**

Door Face Dimensions: 15" wide, 13 1/2" high, 1" thick.

**Sliding Tray**

Dimensions: 12" wide, 2" high, 25 1/2" deep.  
 Weight Capacity: 75 lbs.

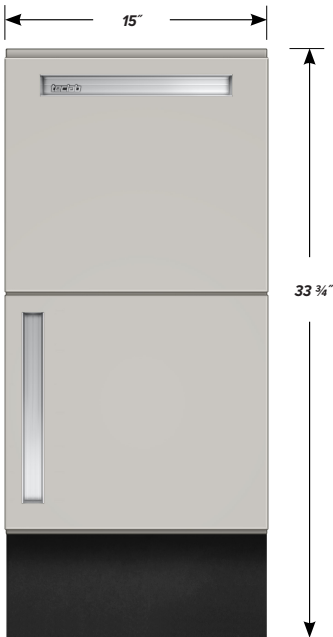

**FRONT**



**SIDE**