

MIXING FAUCET WITH VACUUM BREAKER FM6-412-VB

BASE FAUCET:

Single hole deck mount mixing faucet with polished chrome plated brass body, polished chrome plated swivel/rigid vacuum breaker nozzle with serrated tip and 0.5 GPM flow control, compression cartridges with spring checks, wrist action handles and 18" flexible stainless steel supply hoses with 1/2" NPSM / 5/8" compression connections. Certified to ASME A112.18.1 / CSA B125.1, NSF 61-Section 9, NSF 372 and ASSE 1001. Meets ADA ANSI/ICC A117.1 requirements.

Body: Single Hole

Body Mount: Deck

Cartridge-Spindle: 1/4 Turn Compression

Feature: Atmospheric VB

Feature: Spring Check

Handle: 4 in Wrist

Material: Brass Construction

Inlet: 5/8 in Copper Tubing or 1/2 in NPSM Male

Nozzle: 6 in Gooseneck Swivel

Nozzle Flow: 0.5 GPM (1.89 L/min)

Nozzle Outlet Type: Flow Regulator

Nozzle Outlet Type: Serrated Tip

Pressure Max: 125 psi (8.62 bar)

Pressure Min: 20 psi (1.38 bar)

Temperature Max: 140 F (60°C)

Temperature Min: 40 F (4.44°C)



SERRATED HOSE END

3/8-18 NPT
Male Thread

2 3/16"
[56mm]


10 Serrations, Accepts
Plastic or Rubber Tubing
Up to 1/2" I.D.

le Thi

Washer



MIXING FAUCET WITH VACUUM BREAKER FM6-412-VB

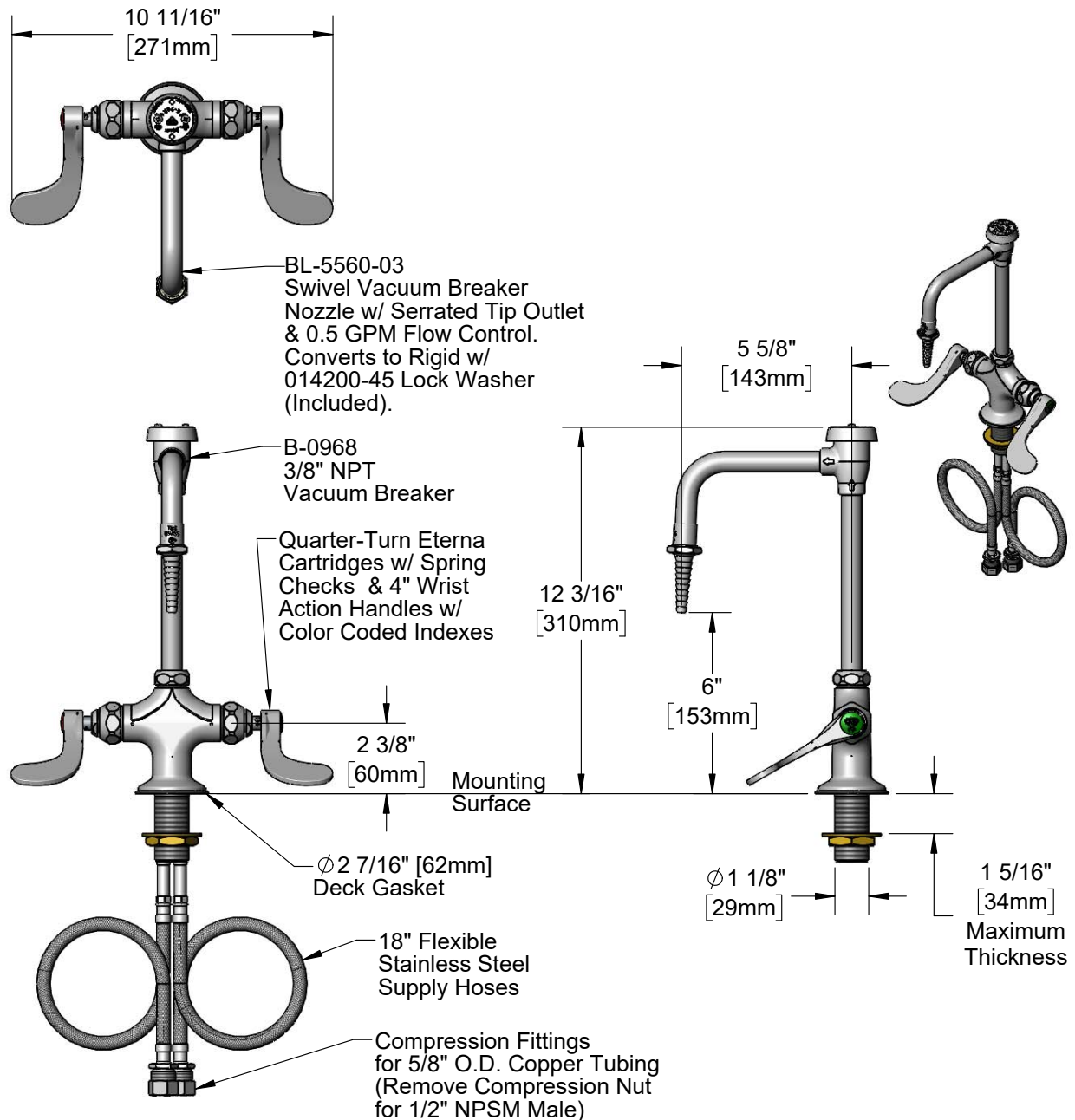

Rough-In Requirement:
Ø 1 1/2" [38mm] Mounting Hole


Product Specifications:

Single Hole Deck Mount Mixing Faucet, Quarter-Turn Eterna Cartridges w/ Spring Checks, 4" Wrist Action Handles, Swivel/Rigid Vacuum Breaker Nozzle, Serrated Tip Outlet w/ 0.5 GPM Flow Control & 18" Flexible Stainless Steel Supply Hoses w/ 1/2" NPSM Male Inlets

Product Compliance:

ASME A112.18.1 / CSA B125.1
NSF 61 - Section 9
NSF 372 (Low Lead Content)
ANSI A117.1 (ADA)
ASSE 1001 (VB)

MIXING FAUCET WITH VACUUM BREAKER FM6-412-VB

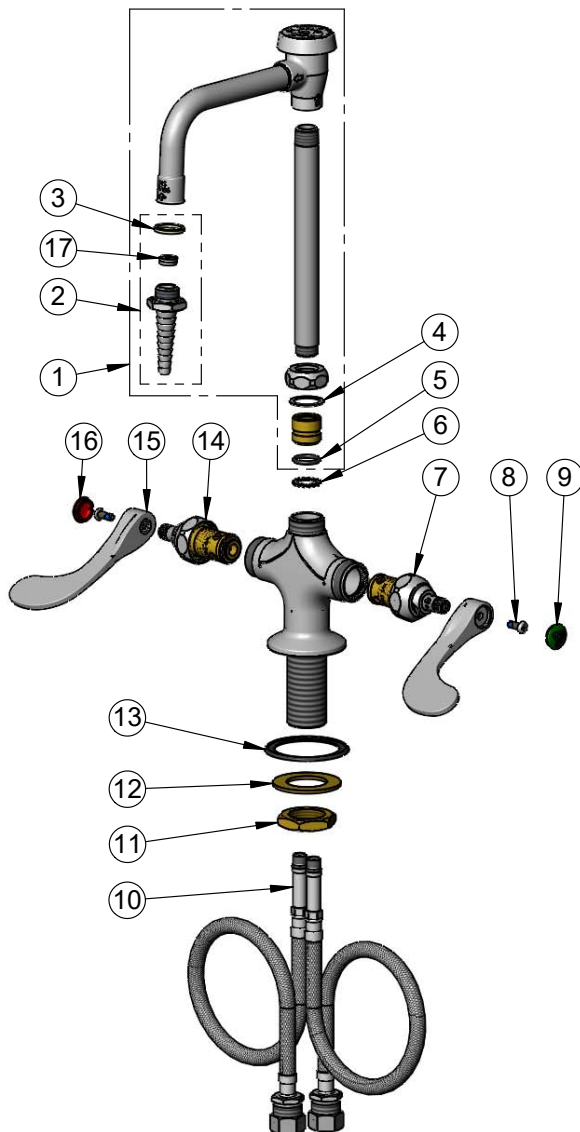

ITEM NO.	SALES NO.	DESCRIPTION
1	BL-5560-03	Swivel Vacuum Breaker Nozzle w/ Serrated Tip & 0.5 GPM Flow Disc
2	B-0198-FD05	Serrated Tip w/ 0.5 GPM Flow Disc & Washer
3	001051-45	Washer
4	009538-45	Swivel Washer
5	001074-45	O-ring
6	014200-45	Star Washer, Anti-Rotation
7	012442-40NS	Quarter-Turn Eterna Cartridge w/ Spring Check, LTC
8	000925-45	Lab Handle Screw
9	001191-45NS	Green "CW" Lab Button Index
10	012534-45	18" Flexible Supply Hose (2)
11	000965-45	Lock Nut
12	002290-45	Lock Washer
13	145L	Black Plastic Gasket
14	012443-40NS	Quarter-Turn Eterna Cartridge w/ Spring Check, RTC
15	B-WH4-NS	Wrist Action Handle - New Style
16	001194-45NS	Red "HW" Lab Button Index
17	FD05	Flow Disc, 0.5 GPM [1.9 L/min.] (Lime)